# Benutzerhandbuch User manual

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DE EN

## AFH Parffino® Paraffinbad

Model FHC-4000D

AFH Webshop®

EN

### 2. Heating Modes

#### A. Manual Heating Mode

- » In manual heating mode, you manually set the desired temperature using the plus and minus buttons (+) / (-), and the device maintains that temperature until switched off.
- Press the plus and minus buttons (+) / (-) to set the target temperature (30–70 °C / 86–158 °F).
- Maintains the set temperature until powered off.

#### B. Smart Heating Mode

- » In smart heating mode, the appliance automatically heats up to the ideal application temperature of 53 °C / 127 °F and keeps it constant without the need for manual settings.

- Tap the Mode Button (Ⓜ) to activate „HEAT“ status.

- Automatically heats to 53 °C / 127 °F; 3 beeps signal readiness.

#### C. Reservation Heating Mode

- » In reservation heating mode, the device starts heating at a delayed time so the wax is ready precisely when you need it.
- With red indicator on, hold the Mode Button (Ⓜ) for 3 seconds to enter reservation.

- Press the plus and minus buttons (+) / (-) to set a countdown timer (2–99 hours).

- Heating activates automatically during the final 2 hours of the countdown.

### 3. Paraffin Therapy Process

- Temperature Test: Use a spatula to dab a small amount of wax onto the inner side of your wrist to test comfort.

#### Immersion:

- Dip hands or feet repeatedly (3 seconds per dip, 5–8 cycles).
- Retain the wax coating after the final dip.

#### Heat Retention:

- Wrap treated areas with plastic film and thermal mitts or boots for 15 minutes.

#### Post-Treatment:

- Peel wax from extremities toward the torso.
- Reuse residual wax up to 3 times, filtering out any debris beforehand.

#### Shutdown & Maintenance

- Press the Power Button (Ⓜ) to turn off the device.
- Unplug the unit and allow it to cool completely.

#### CLEANING AND MAINTENANCE

- Turn off the unit and allow the wax to cool completely.
- Remove any remaining wax residue with the silicone spatula.
- Wipe the exterior of the unit with a dry cloth. Never immerse the device in water!

#### SAFETY AND WARNINGS

1. If the power cord is damaged, it must be replaced by the manufacturer, its service agent, or a similarly qualified person to avoid hazards.
2. This symbol indicates that this product must not be disposed of with regular household waste in the EU. To avoid environmental or health damage from uncontrolled waste disposal, recycle this product responsibly and support the sustainable reuse of materials. Please use appropriate return and collection systems or contact the retailer where the product was purchased for proper disposal.
3. This device may be used by children aged 8 and above, and by persons with reduced physical,

EN

sensory or mental capabilities, or lack of experience and knowledge, if they are supervised or instructed in the safe use of the device and understand the associated hazards.

4. Children must not play with the device.

5. Cleaning and maintenance must not be carried out by children without supervision.

6. Note: When using electrical devices – especially around children – basic safety precautions should always be followed. Read all instructions carefully before using the wax bath.

7. This device is intended for heating paraffin wax only.

8. To avoid burns, electric shocks, fire, or injury, do not use the device if:
  - The power cord or plug is damaged
  - The device itself is damaged
  - The device is not functioning properly

- The power cord or plug is damaged
- The device itself is damaged
- The device is not functioning properly

9. Keep the cord away from hot surfaces. Always operate the device on a flat surface. Do not immerse the device in water or other liquids. Keep the area around the unit dry. Do not place the device where it could fall or be pulled into a bathtub or sink. If the unit falls into water, do not touch it! Unplug it immediately. Do not wrap the cord around the unit during operation. Unplug the unit after use. Allow it to cool completely before handling, cleaning, or storing. Keep out of reach of children. The device remains hot for some time after use.

#### WARRANTY CONDITIONS

The AFH Parffino® Paraffin Bath XL Model FHC-4000D comes with a 1-year warranty from the date of purchase. This guarantee applies to all faulty parts that are no longer functional due to material or manufacturing defects. The proof of purchase serves as proof and is required to claim under the guarantee.

The AFH Webshop wishes you a pleasant and relaxing experience with your paraffin wax bath!

**AFH Parffino® Paraffinbad**  
Model: FHC-4000D | Power: 200W  
Voltage: 220–240 V 50/60Hz

Guangzhou Fourto Sanitary Products Co., Ltd. No.7, Shiyuan Rd 2, Xinhua Street, Huadu District, Guangzhou, China.

**AFH-Webshop®** · Inhaber Rainer Zumhasch  
Pyrmonter Straße 50 · 32676 Lügde / NRW · Germany


**AFH Webshop®** · Inhaber Rainer Zumhasch  
Pyrmonter Straße 50 · 32676 Lügde / NRW  
· Germany

**AFH Parffino® Paraffinbad**  
Model: FHC-4000D | Power: 200W  
Voltage: 220–240 V 50/60Hz

Guangzhou Fourto Sanitary Products Co., Ltd. No.7, Shiyuan Rd 2, Xinhua Street, Huadu District, Guangzhou, China.

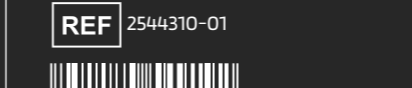**AFH-Webshop®** · Inhaber Rainer Zumhasch  
Pyrmonter Straße 50 · 32676 Lügde / NRW · Germany

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**EN USER MANUAL**

This user manual includes **product specifications, safety warnings, operating instructions, cleaning and maintenance tips**, as well as important **safety and caution notices**.

Product Name:	AFH Parffino® Paraffin Bath
Model Number:	FHC-4000D
Power Rating:	200 W
Voltage Compatibility:	US Standard: 110-120V~ 60Hz / EU Standard: 220-240V~ 50Hz
Temperature Range:	30-75°C (86-167°F)
Capacity:	4500 ml
Dimensions:	408 mm (L) x 254 mm (W) x 170 mm (H)
Construction:	Body: Food-Grade ABS Resin, Lid: Heat-Resistant PS Material
Certifications:	CE (EMC/LVD), FCC, RoHS, Reach
Electrical Class:	Class I

**SAFETY WARNINGS**

Please read all instructions carefully before first use.

**ELECTRICAL SAFETY**

- Use the device only with the specified voltage.
- Do not operate the unit if the power cord or plug is damaged. Contact qualified service personnel for repairs.
- Keep the power cord away from hot surfaces and liquids.

**USER SAFETY**

- Risk of Burns: Never touch the wax or insert body parts into the device during the heating process.
- Children and Vulnerable Users:

- » Supervise children aged 8 years and older, and persons with reduced physical, sensory, or mental capabilities.

- » Keep the device out of reach of children, especially when it is hot.


**ENVIRONMENT:**

- Place the device on a flat, dry surface. Never immerse in water.
- If the device falls into liquid, unplug it immediately.

**OPERATIONAL SAFETY**

- This device is intended for heating paraffin wax only.
- Allow the unit to cool completely before cleaning or storage.
- Dispose of the unit according to electronic waste guidelines (EU WEEE Directive).

**OPERATING INSTRUCTIONS**

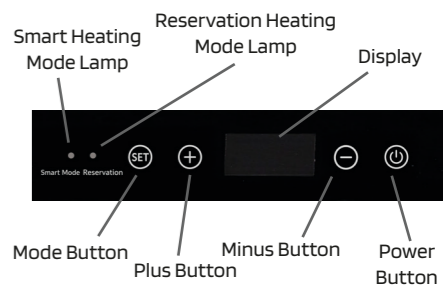**1. Setup and Initial Use**

- Place the warmer on a stable, dry surface.
- Connect to power and press the Power Button (ON/OFF) to activate.
- Open the lid and add cosmetic-grade paraffin wax (≤4 L recommended).