Benutzerhandbuch User manual

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DE EN


AFH Parffino® Paraffinbad XL

Model FHC-9000A

AFH Webshop®

EN OPERATING INSTRUCTIONS

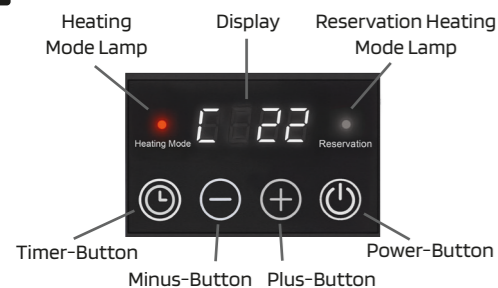1. Setup and Initial Use

- Place the device on a stable, dry surface.
- Remove the lid and fill the silicone container with paraffin wax.
- Secure the lid back in place.
- Connect the device to power and press the Power Button to switch it on.

2. Heating Modes

A. Manual Heating Mode

» In manual heating mode, you manually set the desired temperature using the plus and minus buttons / , and the device maintains that temperature until switched off.

Set the desired temperature (30 to 75 °C or 86°F to 167°F) using the plus and minus buttons / .

Heating mode lamp:
A red light indicates that the heating is active, while a green light signals that the set temperature has been reached.

The device maintains the set temperature until powered off.

B. Reservation Heating Mode

» In reservation heating mode, the device starts heating at a delayed time so the wax is ready precisely when you need it.

Tap the Timer Button until the Reservation

Heating Mode light turns red to activate (default: H04 = 4 hours).

Set the desired duration using the plus and minus buttons / (4–99 hours).

Heating starts 4 hours before the set time ends.

Completion:

- Three audible beeps signal that the wax is ready.
- The device automatically cools to 53 °C (127°F) and maintains this temperature.

3. Paraffin Therapy Process

- Preparation: After the wax has completely melted, insert the heat-insulating grid.
- Temperature Test: Briefly dip fingertips or toes to test the wax temperature.

Application:

Immerse hand or foot for 3 seconds → Remove for 3 seconds → Repeat 5–8 times.

Apply a final layer of wax with a brush.

Heat Retention:

Cover the treated area with a plastic liner and thermal mitt for about 15 minutes.

Post-Treatment:

Gently peel off the hardened wax from the wrist or ankle.

Dispose of the used wax and liners.

4. Shutdown & Maintenance

Press the Power Button to turn off the device.

Unplug the unit and allow it to cool completely.

EN CLEANING AND MAINTENANCE

Clean the silicone tray before first use.

Thanks to the „Easy Clean Pro“ function, cleaning the wax bath is particularly easy:

Turn off the unit and allow the wax to cool completely.

Remove the silicone container along with the wax and empty it. Solidified wax can be reused.

Remove any remaining wax residue with the silicone spatula.

Wipe the interior with medical alcohol or a silicone cleaner and a soft cloth.

Wipe the exterior of the unit with a dry cloth. Never immerse the device in water!

Tip:

By using multiple silicone containers, you can easily switch between different wax aromas or scents without needing to clean the device each time.

SAFETY AND WARNINGS

1. If the power cord is damaged, it must be replaced by the manufacturer, its service agent, or a similarly qualified person to avoid hazards.

2. This symbol indicates that this product must not be disposed of with regular household waste in the EU. To avoid environmental or health damage from uncontrolled waste disposal, recycle this product responsibly and support the sustainable reuse of materials. Please use appropriate return and collection systems or contact the retailer where the product was purchased for proper disposal.

3. This device may be used by children aged 8 and above, and by persons with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, if they are supervised or instructed in the safe use of the device and understand the associated hazards.

Children must not play with the device.

Cleaning and maintenance must not be carried out by children without supervision.

Note: When using electrical devices – especially around children – basic safety precautions should always be followed. Read all instructions carefully before using the wax bath.

This device is intended for heating paraffin wax only.

To avoid burns, electric shocks, fire, or injury, do not use the device if:

- The power cord or plug is damaged
- The device itself is damaged
- The device is not functioning properly

9. Keep the cord away from hot surfaces. Always operate the device on a flat surface. Do not immerse the device in water or other liquids. Keep the area around the unit dry. Do not place the device where it could fall or be pulled into a bathtub or sink. If the unit falls into water, do not touch it! Unplug it immediately. Do not wrap the cord around the unit during operation. Unplug the unit after use. Allow it to cool completely before handling, cleaning, or storing. Keep out of reach of children. The device remains hot for some time after use.

WARRANTY CONDITIONS

The AFH Parffino® Paraffin Bath XL Model FHC-9000A comes with a 1-year warranty from the date of purchase. This guarantee applies to all faulty parts that are no longer functional due to material or manufacturing defects. The proof of purchase serves as proof and is required to claim under the guarantee.

The AFH Webshop wishes you a pleasant and relaxing experience with your paraffin wax bath!

AFH Webshop® · Inhaber Rainer Zumhasch
Pyrmonter Straße 50 · 32676 Lügde / NRW
· Germany

AFH Parffino® Paraffinbad XL
Model: FHC-4000D | Power: 350W
Voltage: 220-240 V 50/60Hz

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USER MANUAL

This user manual includes **product specifications, safety warnings, operating instructions, cleaning and maintenance tips**, as well as important **safety and caution notices**.

Product Name:	AFH Parffino® Paraffin Bath XL
Model Number:	FHC-9000A
Power Rating:	350 W
Voltage Compatibility:	US Standard: 110-120V AC / EU Standard: 220-240V AC
Temperature Range:	30-75°C (86-167°F)
Capacity:	9000 ml
Dimensions:	420 mm (L) x 300 mm (W) x 200 mm (H)
Construction:	201 Stainless Steel Tank, ABS Frame, PC Lid
Certifications:	CE (EMC/LVD), FCC, RoHS, Reach
Electrical Class:	Class I

SAFETY WARNINGS

Please read all instructions carefully before first use.

ELECTRICAL SAFETY

- Use the device only with the specified voltage.
- Do not operate the unit if the power cord or plug is damaged. Contact qualified service personnel for repairs.
- Keep the power cord away from hot surfaces and liquids.

USER SAFETY

- Risk of Burns: Never touch the wax or insert body parts into the device during the heating process.
- Children and Vulnerable Users:
 - » Supervise children aged 8 years and older, and persons with reduced physical, sensory, or mental capabilities.
 - » Keep the device out of reach of children, especially when it is hot.

ENVIRONMENT:

- Place the device on a flat, dry surface. Never immerse in water.

- If the device falls into liquid, unplug it immediately.

OPERATIONAL SAFETY

- This device is intended for heating paraffin wax only.
- Allow the unit to cool completely before cleaning or storage.
- Dispose of the unit according to electronic waste guidelines (EU WEEE Directive).