MASSAGELIEGE

Massage table • Table de massage • Lettino per massaggi
Masažni ležalnik • Stol za masažu

- Bedienungsanleitung
- Instruction manual
- Mode d'emploi
- Istruzioni per l'uso
- Navodila za uporabo
- Upute za uporabu

Deutsch 2
English 16
Français 30
Italiano 44
Slovenščina 58
Hrvatski 72

RB#65620 (Typ IB65584)
(CREME)

2017-07

**COMFORT
MED**

Royalbeach®
SPORTS & ACTION

Inhalt

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Inhalt




Overview.....	18
Scope of delivery /equipment parts	18
General	20
Reading and keeping instruction manual	20
Legend.....	20
Safety	21
Intended use	21
Safety instructions	21
Checking massage couch and scope of delivery	22
Use	22
Assembling massage couch	22
Setting lying height	23
Attaching armrests	24
Attaching head rest and hand rest	24
Lying down on massage couch	25
Standing up from massage couch	25
Folding massage couch	26
Cleaning	27
Maintenance	27
Aufbewahrung	27
Technical data	28
Disposal	28
Disposing of packaging	28
Disposing of massage couch	28
Hotline.....	28









Scope of delivery / Overview

- 1 Lying surface
- 2 Plastic foot, 4 × at both sides
- 3 Carrying handle
- 4 Support leg, 4×
- 5 Adjusting screw, two at each support leg
- 6 Wire cable
- 7 Locking, 2×
- 8 Armrest, 2×
- 9 Head rest pad, removable
- 10 Hand rest locking
- 11 Pad of the head rest
- 12 Locking lever of the head rest
- 13 Frame of the head rest
- 14 Carrier bag
- 15 Closure, 2×
- 16 Rubber mounting (for the head rest and the pad)
- 17 Rubber mounting (for the armrests and the hand rest)

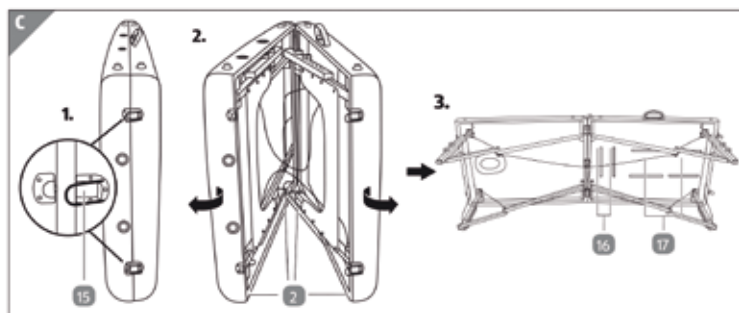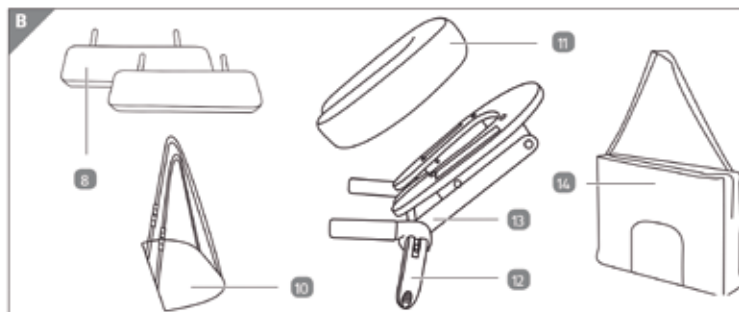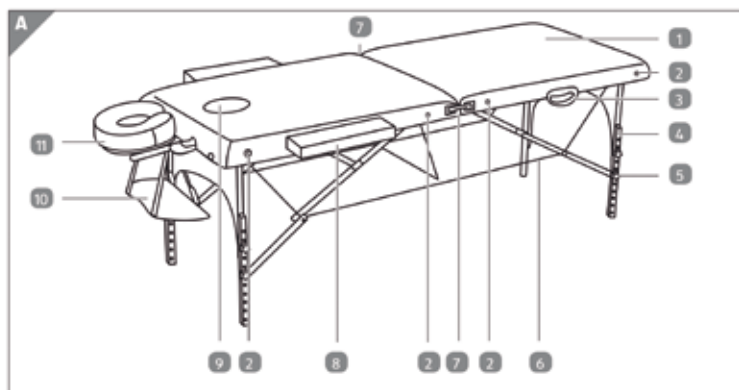Overview

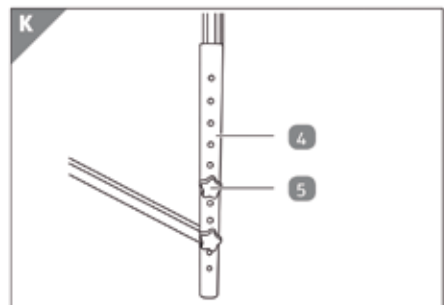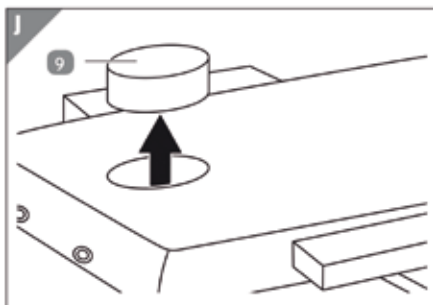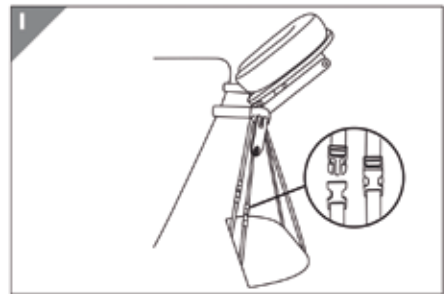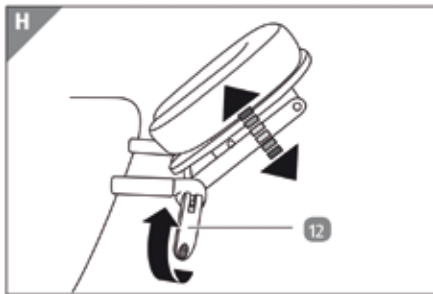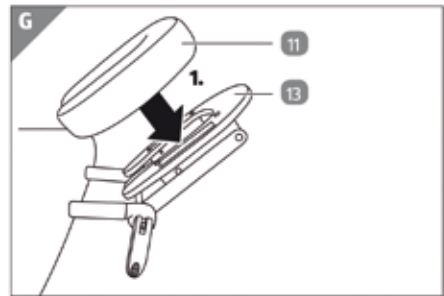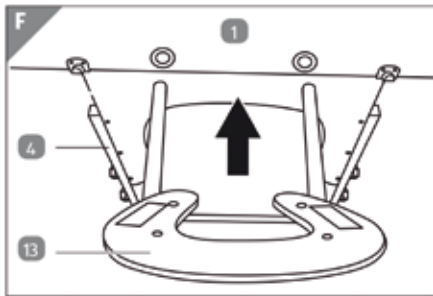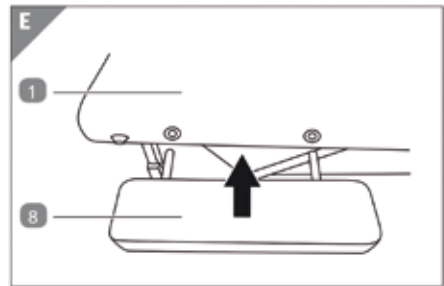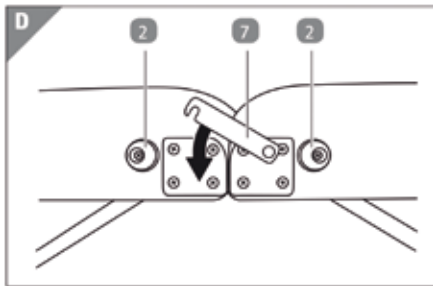General

Reading and keeping instruction manual

This instruction manual belongs to this massage couch. It contains important information on the mounting and the handling. Please read the instruction manual, especially the safety instructions, carefully before using the massage couch. Non-compliance with this instruction manual can lead to severe injuries or damage to the massage couch.

The instruction manual is based on the valid norms and regulations in the European Union. Please also observe country-specific directives and laws abroad. Please keep this instruction manual for further use. If the massage couch is handed on to third parties, please also hand on this instruction manual.


Legend

The following symbols and signal words are used in this instruction manual, on the massage couch or on the packaging.


WARNING!

This signal word refers to a danger with a medium risk level which, unless avoided, can result in death or a severe injury.

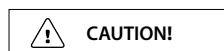CAUTION!

This signal word refers to a danger with a low risk level which, unless avoided, can result in a slight or moderate injury.

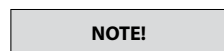NOTE!

This signal word warns of possible material damages.

This symbol provides you with useful additional information on the assembly or the operation.

www.tuv.com
ID 000038655

The TÜV [MOT] seal confirms that the massage couch is safe in the case of intended use. The test seal GS (Geprüfte Sicherheit [Tested for Safety]) confirms the conformity of the massage couch with the Produktsicherheitsgesetz (ProdSG) [German Product Safety Act].

The wood used in the product comes from FSC® certified forest holdings which comply with the FSC® principles. The Forest Stewardship Council® (FSC) aims at a worldwide responsible forest management.

Safety

Intended use

The massage couch is exclusively designed as a lying surface. It is exclusively intended for private use and is not suited for the commercial sector. The massage couch is not a toy.

Please use the massage couch only in accordance with the descriptions in this instruction manual. Any other use is regarded as not intended and can lead to material damages or even personal injury.

The manufacturer or dealer does not take any liability for damages which are the result of not intended or wrong use.

Safety instructions

Please read the safety instructions carefully, before assembling and using the massage couch.

CAUTION!

Risk of injury!

Improper use can lead to risk of injury!

- Assemble the massage couch horizontally and wobble-free on an even, solid ground, so that it cannot wobble and tip over.
- Always secure the massage couch first with the lockings on both sides, before lying on the massage couch.
- Make sure that no children play with the massage couch.
- Carefully fold the massage couch open and shut in order not to crush your fingers. Moveable parts generally bear the danger of crushing.

NOTE!

Risk of damage!

Improper use of the massage couch can lead to damages of the massage couch.

- Do not expose the massage couch to direct sunlight.
- Do not place the massage couch near sources of heat.
- Do not place the massage couch near hot surfaces.
- Load the massage couch with max. 125 kg user weight.
- No longer use the massage couch, if the support legs or the lying surface have cracks or the wire cables fray.
- Replace damaged parts only with suitable original spare parts.
- Have the repairs carried out only by qualified personnel.
- Do not change the couch's construction.
- Do not change the position of the moveable parts.
- Place a blanket under the massage couch when placing it on its side (e.g. during installation or when extending / shortening the support legs). Checking massage couch and scope of delivery

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WARNING!

Danger of suffocation!

If children or animals play with the packaging foil, they can get caught in it and suffocate.

- Do not allow children and animals to play with the packaging foil.

NOTE!

Risk of damage!

If you open the packaging carelessly with a sharp knife or other sharp objects, the massage couch can quickly be damaged.

- Therefore be very careful when opening it.

1. Remove the massage couch from the packaging.
2. Check whether the delivery is complete (see fig. A and B).
3. Check whether the massage couch or the individual parts are damaged.

If this is the case, do not use the massage couch. Please contact the manufacturer via the service hotline indicated on the warranty card and then send back the defective massage couch, if necessary.

Use

The workload of the massage couch results from the body weight of the lying person (user weight) and the maximally exercised pressure during the massage. The massage couch is designed for carrying a person with a total weight of up to 125 kg. The pressure during the massage must not exceed 100 kg. This results in the total maximum workload of 225 kg.

Lie down in the middle of the massage couch in order to evenly distribute the workload.

Assembling massage couch

CAUTION!

Risk of injury!

Improper use can lead to risk of injury!

- Place the massage couch horizontally and wobble-free on an even, solid ground so that it cannot wobble and tip over.
- Always secure the massage couch first with the locking on the side, before lying on the massage couch.
- Carefully fold the massage couch open in order not to crush your fingers.
- Check the massage couch before each use for e.g. loose screws, cracks and other things. Do not use the massage couch if it is damaged.

NOTE!**Risk of damage!**

If you do not have enough moving space when assembling the massage couch you can damage adjacent furniture and objects.

- Assemble the massage couch with a minimum clearance of 2 m on all sides.
 - Do not open the massage couch by force if the support legs are stuck.
Search for and remove the cause of this failure and only then assemble the massage couch afterwards.
1. Unpack the massage couch from the carrier bag **14**.
 2. Position the closed massage couch at a suitable place on the four plastic feet **2** (see step 1 in **fig. C**).
 3. Open both closures **15** and fold the massage couch open (see step 2 in **fig. C**).
 4. Remove all individual parts from the rubber mountings **16** and **17**.
 5. Completely unfold the two parts of the lying surface **1** (see step 3 in **fig. C**).
Thus the support legs unfold automatically in a 90° angle.
 6. Lift the massage couch in the middle and place it on the support feet.
 7. Make sure that the lying surface is even. If the lying surface is not yet even, adjust the wire cables **6** and apply little pressure on the lying surface, if necessary.
 8. Secure the two halves of the lying surface with the lockings **7** (see **fig. D**).

You have assembled the massage couch successfully and now can adjust the lying height and the armrests, as well as attach the head and arm supports.

Setting lying height You can adjust the height of the massage couch by loosening and adjusting the adjustment screws at the support legs.

Adjust the lying height before attaching the armrests, the head rest and the hand rest.


1. Place the massage couch on its left side.
2. Loosen the two adjustment screws **5** at the upper support legs **4**.
3. Remove the adjustment screws and the washers and put them aside.
4. Remove the upper wooden strip with the boreholes upwards from the two screws.
5. Evenly adjust the position of the support legs by placing the wooden strip with the boreholes on the support legs so that both screws are led into one borehole each.
Both support legs must be equally extended or shortened.
6. Attach the washers to the screws.
7. Tighten the adjustment screw clockwise.

You have successfully extended / shortened the support legs of the left side.

8. Place the massage couch on its right side.
9. Repeat step 2 to 8 on the right side. Please make sure that all of the four support legs have the same length at the end

You have successfully extended / shortened the support legs of the massage couch.

Attaching armrests

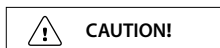Risk of injury!

If you lean on the armrests, these can fall off because they are not firmly fixed onto the massage couch. You can injure yourself.

- Do not lean on the armrests, but only use them to rest your arms on.
- Push the armrests 8 all the way into the lateral holes of the lying surface **1** (see **fig. E**).

Attaching head rest and hand rest

Risk of injury!

If the locking lever of the head rest is not fixed correctly, the head rest can suddenly tilt down. This can lead to injury of the person on the massage couch.

- Correctly fix the locking lever of the head rest before using the massage couch. The head rest consists of a frame and a pad.

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This is how to attach the head rest:

1. Insert the frame of the head rest **13** in the two holes at the head end of the lying surface **1** (see **fig. F**).
2. Fix the frame of the head rest at the desired angle by turning the locking lever of the head rest **12** (see **fig. H**).
3. Attach the pad of the head rest **11** by means of the velcro fastener to the frame (see **fig. G**).

You have successfully attached the head rest to the massage couch and can now attach the hand rest locking **10**.

This is how to attach the hand rest locking to the head rest:

- Hook the hand rest locking into the frame of the head rest (see **fig. I**). Now you can set the desired length and position of the hand rest locking.

Lying down on massage couch

CAUTION!

Risk of injury!

If you lie down incorrectly on the massage couch, it will be loaded incorrectly. This can result in the massage couch tipping over and you could injure yourself.

- Do not sit down on the head end or foot end of the massage couch.
- Always position yourself in the middle of the lying surface.

This is how to lie down correctly on the massage couch:

1. Sit down on the middle of the lying surface 1.
2. Lie down in the desired position on the massage couch. Make sure that you always position yourself in the middle of the massage couch in order to evenly spread the working load.

Standing up from massage couch

CAUTION!

Risk of injury!

If you stand up incorrectly from the massage couch, it will be loaded incorrectly. This can result in the massage couch tipping over and you could injure yourself.

- Do not sit on the head end or foot end of the massage couch when standing up.

This is how to safely stand up from the massage couch:

1. Sit up in such a way that you are sitting in the middle of the lying surface 1.
2. Stand up slowly.

Quick movements are not recommended after longer relaxed lying, because your blood circulation has to adjust to this first.

Folding massage couch

CAUTION!

Risk of injury!

Improper use can lead to risk of injury!

- Carefully place the massage couch on the side and fold it with caution in order not to crush your fingers at the moveable parts.

NOTE!

Risk of damage!

If you do not have enough moving space when folding the massage couch, you can damage adjacent furniture and objects.

- Fold the massage couch with a minimum clearance of 2 metres on all sides.
- Only fold the massage couch together if it is clean and dry.

If you have extended the length of the support legs, you might have to shorten them again before folding away the massage couch.

1. Remove the armrests **8**, the head rests **11 / 13** and the hand rest **10** from the massage couch.
2. Loosen the lockings **7** at both sides of the lying surface **1** (see step 3 in **fig. C**).
3. Place the massage couch on the left side. The massage couch is now standing on the four plastic feet **2**.
4. Store the individual parts with the help of the rubber mounts **16 / 17** under the massage couch.
5. Slowly fold the support legs back under the massage couch. Make sure that the wire cables **15** do not get tangled up and that the individual parts do not impede the support legs.
6. Fold the two parts of the lying surface together and close the closings **15**. Make sure that the massage couch can be easily folded. It must not be under tension.

You have successfully folded the massage couch.

Cleaning

NOTE!**Risk of damage!**

If you use unsuitable cleaning agents or utensils, these can damage the surfaces of the massage couch and the individual parts.

- Do not use aggressive, acid-, vinegar or bleach-containing cleaning agents, impregnating agents or other chemicals.
- Do not use brushes with metal or nylon bristles.
- Do not use any sharp or metallic cleaning accessories such as knives, hard palette knives and the like.

This is how to remove slight contaminations:

1. Wipe the lying surface 1 with a slightly moistened cloth and some mild soap, if necessary.
2. Afterwards dry all parts completely with a dry, lint-free cloth, before re-using the massage couch or folding it away and storing it in the bag

This is how to remove stronger contaminations:

1. Wipe off stronger contaminations with a damp cloth and some mild detergent.
2. Wipe off the affected area with warm, clear water.
3. Afterwards dry all parts completely with a dry, lint-free cloth, before re-using the massage couch or folding it away and storing it in the bag. We recommend to have stubborn contaminations cleaned professionally.

We recommend to have stubborn contaminations cleaned professionally.

Maintenance

**CAUTION!****Risk of injury!**

The use can lead to the gradual loosening of screw connections. If you do not fasten these in time, the massage couch can fall apart and you can injure yourself.

- Carry out regular visual inspections with regard to cracks and loose screws. Do not use the massage couch before having tightened loose screws and having had defective parts replaced by a specialist.
- Check on a regular basis whether screw connections have become loose.
- Tighten loose screws.

Storage

- Store the massage couch in folded position in the carrying bag **14**
- Store the massage couch in a cool, dry and clean place.

Technical data

Model:	RB#65620 (Typ IB65584) (creme)
Assembly (L × W × H):	2170 × 920 × 635 - 860 mm
Dimensions of the lying surface (L x W):	1,910 × 700 mm
Weight including accessories:	approx. 15.2 kg
Max. user weight:	125 kg
Max. work load:	225 kg
Material of the frame:	solid wood beech
Bracing:	wooden cross beams and wire cables
Material of the cover:	oil- and water-resistant PU
Height adjustment:	from 63.5 to 86.0 cm

Disposal

Disposing of packaging

Dispose of the packaging with correct sorting. Place cardboard and paperboard in the waste paper, foils in the collection of reusable materials.

Disposing of massage couch

Dispose of the massage couch in accordance with the provisions and laws which are valid in your country and / or your region.

Hotline

For technical inquiries regarding this equipment, please use our service hotline or send us an email.

- Telephone: +49/ 86 85/98 89 80 1
- Mail: service@royalbeach.de
- Spare part order: Indicate exact name and order number!
- Name of equipment:: Massage table
- Item no.: RB#65620 (Typ IB65584) (creme)

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